

PARAMATER

Material: Silk
Heating System: Graphene Heating System
Rated Power: 5W
Rated Voltage: 5V
Output Current: 0.95A
Temperature: 40°C - 60°C (104°F - 140°F)

APPLICATION

The heating system can effectively activate cells, improve blood circulation, promote metabolism and balance secretion.
Perfect to reduce pain in the neck and cervical vertebrae, and relieve symptoms such as neuralgia and muscle soreness.

WASHING

- 1. The ARRIS heating neck wrap supports hand wash.
- 2. To protect the heating system, please wash the neck wrap at 30°C and at a gentle spin.



HOW TO USE

1. Insert the DC plug into the port on neck wrap.

2. Connect the USB plug to any 5V power supply. (It can supply power from 5V USB, such as power bank, PC, car USB adapter, etc.)
(5V Power Bank Not Included.)

3. Press ON/OFF button shortly on the controller to turn on/off the heating system.


4. Wear heating neck wrap, hot compress time is recommended within 30 minutes.

3.1 To adjust the temperature setting
Press the HOT button shortly on the controller to choose the desired temperature.

3.2 To adjust the timer setting
Press the TIME button shortly on the controller to choose the desired time.

TIME AND TEMPERATURE SETTING

There are three adjustable temperature/ three time settings, the temperature and time are quite distinct. Easy to use.

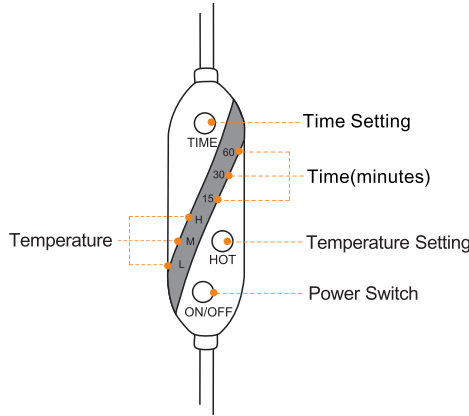

Temperature:

- H: High 60°C (140°F)
- M: Medium 50°C (122°F)
- L: Low 40°C (104°F)

NOTE: The temperature is directly tested from the internal heating element, that means the external temperature you feel when wearing it will be slightly lower.

Time:

- 15 minutes
- 30 minutes
- 60 minutes


NOTE: Recommended time: 30 minutes

! PRECAUTION

- 1. DO NOT wear the heating neck wrap when go into deep sleep at night. The preset longest heating time is 60 minutes. If exceed this time frame, customers need to restart the neck wrap again. Suggest time is 30 minutes.
- 2. Due to the difference in personal skin sensitivity, very few users can be allergic to temperature. If the users fell unwell, please discontinue to use temporarily to avoid burns or skin allergies.
- 3. Never use the neck wrap if it is wet or moist.
- 4. Do not put explosives, combustible and combustible objects near the heating element.
- 5. It takes around 3-5 minutes to heat up. It feels warm, not hot. Too hot might cause hurt to your skin.
- 6. Many 5V 1A mobile batteries in the market only have a current output of only 600-700mah, so they cannot provide enough power for the rear cover.Please use 5v 2a power supply when the temperature is too low.
- 7. Those with pacemakers, pregnant women, young children and those who can not express their willingness, please follow the doctor's instructions to us this product.

MULTIPLE WAYS TO POWER THE NECK WRAP

Any 5V/1-2A USB outlet can supply power to our neck wrap.
5V is a safe voltage, no need to worry about the electric leakage problem.

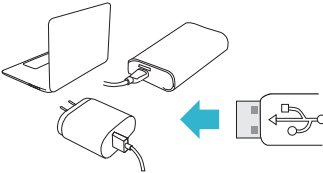

30 DAY MONEY BACK GUARANTEE

We know you will be impressed with our product and customer service, but if you have any questions or concerns, please contact us immediately at support@arrislife.com.



ARRIS
ELECTRIC HEATED NECK WRAP
USER MANUAL

